


Lasermate Group, Inc.

The Friend of Lasers

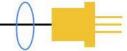

DPF355 Series 355nm UV CW DPSS Laser System up to 10mW

Overview

The DPF355 series diode pumped solid state (DPSS) laser is ideal for applications that require a wavelength of 355nm and output power levels up to 10mW. The laser features a compact design, long lifetime, easy operation, and FDA-compliant system with driver. The laser is widely used in DNA sequencing, flow cytometry, optical instrument, spectrum analysis, measurement, cell sorting, physics experiment, laser rapid modeling, and many other applications.



Specifications

Model Number		DPF355-XYZ
Wavelength (nm)		355±1
Output power (mW)		>3 (X=3), >5 (X=5), >10 (X=10)
Transverse mode		Near TEM ₀₀
Operating mode		CW
Power stability (rms, over 4 hours)		<10% (Z=B)
Beam diameter at aperture (mm)		<2.0
Beam divergence, full angle (mrad)		<1.5
M ² factor		<2.0
Polarization ratio		>50:1
Warm-up time (minutes)		<10
Pointing stability after warm-up (mrad)		<0.05
Beam height from base plate (mm)		45
Spectral purity		>99%
Operating temperature (°C)		10~35
Dimensions of laser head (mm)		211(L)×88(W) ×74(H) mm ³
Weight of laser head (kg)		2.0 kg
Power supply options		
High Power Elite Power Supply (Y=H) :	Input voltage	90-264VAC
	Dimensions	236(L) ×145(W) ×104(H) mm ³
	Weight	2.3 kg
	Notes	CW mode with complete FDA compliant features, such as the turnkey switch and interlock, and easy to operate
High Power Laboratory Power Supply (Y=M)	Input voltage	90-264VAC
	Dimensions	277(L) ×145(W) ×106(H) mm ³
	Weight	2.6 kg
	Notes	Complete FDA compliant features, such as the turnkey switch and interlock, and many other functions; CW mode with adjustable output power knob and operating current LED display
Expected lifetime (hours)		5,000
Warranty period		10 months
FDA Compliance		FDA CDRH Title 21 CFR 1040.10/11 Class IIIb

Remarks:


• Specifications of the CW laser is based on the laser performance at full power output after the specified warmup period. The stability of output power may change when output power is adjusted at a different power level.

Tel: (909)718-0999 | Fax: (909)718-0998 | E-mail: info@lasermate.com | URL: http://www.lasermate.com




Lasermate Group, Inc.

The Friend of Lasers

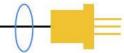

Outline Dimensions (unit: mm)

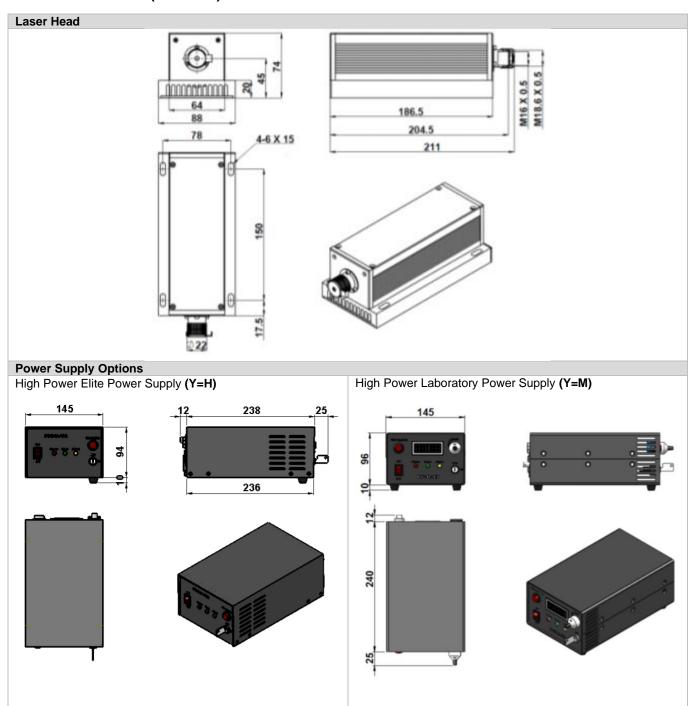

Note: The above specifications are subject to change without notice.

Tel: (909)718-0999 | Fax: (909)718-0998 | E-mail: info@lasermate.com | URL: http://www.lasermate.com